

THE ORGANOID HANDBOOK

BUILDING BETTER ORGANOIDS



biotechne[®]

INTRODUCTION

An organoid is a miniaturized version of an organ produced *in vitro* that shows realistic micro-anatomy, is capable of self-renewal and self-organization, and exhibits similar functionality as the tissue of origin. Organoids are model systems that, in conjunction with advances in cell reprogramming technology and gene editing methods, allow unprecedented insight into human development, disease modeling, drug screening, and transplantation.

Organoids can be classified into those that are tissue-derived and those that are pluripotent stem cell-derived. Tissue-derived organoids typically originate from adult tissues while stem cell-derived organoids are established from embryonic (ESC) or induced pluripotent stem cells (iPSC). Researchers have devised methods to generate physiologically relevant organoid models for many organs, including the intestines, lung, brain, liver, lung, pancreas, and heart. While methods for generating organoids are still evolving, presently they are providing exciting and more accurate systems that are advancing our understanding of basic organ biology and tissue regeneration.

This handbook provides a resource for key publications, protocols, reagents, and troubleshooting recommendations for organoid cell culture.

TABLE OF CONTENTS

Introduction	2
Organoid Culture	3
Tissue-specific Organoids	4
Intestinal.....	4
Gastric	5
Liver	6
Lung.....	7
Brain.....	8
Kidney.....	9
Heart.....	10
Mammary	10
Pancreatic.....	11
Inner Ear.....	12
Organoid Harvesting	13
Organoid Cryopreservation	13
Imaging Organoids	14
Single Cell <i>In Situ</i> Hybridization in Organoids.....	14
Organoid Viability	14

ORGANOID CULTURE

While different methods, such as the use of low adhesion round bottom dishes and bioreactors, have been employed for organoid generation, generally organoids are cultured in tissue culture plates while embedded in “dome” of purified extracellular matrix hydrogels and submerged in organoid-specific culture medium. Multiple organoids are often cultured in one “dome” and, with media changes, submerged organoids can remain in long-term culture to accommodate developmental and maturation timelines.

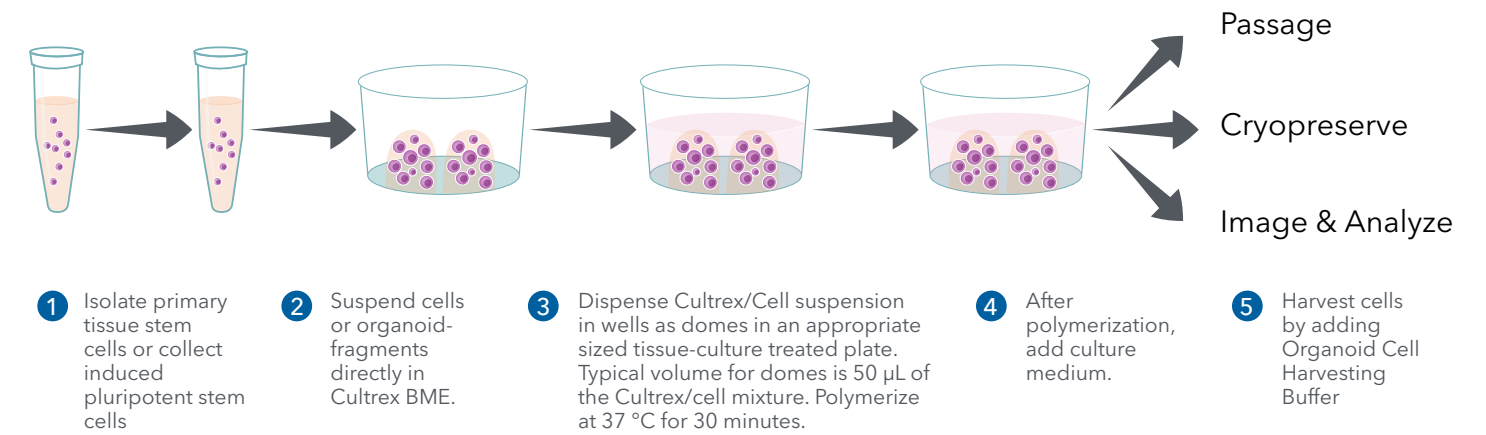


FIGURE 1: General Schematic of Organoid Cell Culture. Individual cells or organoid fragments are embedded within a liquid extracellular matrix (ECM) and dispensed as small droplets onto the surface of a warm tissue culture plastic vessel. The ECM will solidify into a gel after incubation at 37°C and can then be covered with culture medium. Organoids will develop within the dome as 3D structures that can be harvested followed by passaging, cryopreservation, or analysis.

With the increased complexity of organoids and their culture and protocols, the risk of aberrant differentiation and culture variability can also increase. For example, the starting materials (iPSCs or adult stem cells/tissue) interaction with the extracellular matrix is critical for physiological development of the organoid. It is important that individual organoids do not come into contact with one another or with cell culture plastic as this can disrupt or advance organoid development. In addition, the quality and consistency of reagents (e.g., recombinant proteins, small molecules, and extracellular matrix hydrogels) are key elements for developing and maintaining a robust and consistent organoid culture protocols.

Leading organoid researchers agree that reproducibility and culture longevity are the biggest challenges facing organoid biology. Some techniques recommended that can improve long-term culture quality and increase model consistency and are: using reproducible reagents, following consistent culturing protocols (media formulations, splitting protocols, and timing), filtering organoids during passaging to facilitate consistent organoid size, and optimizing organoid density in matrices while ensuring tissue remains fully embedded and does not contact cell culture plastic.


Learn more in our Virtual Organoid Symposium | bio-techne.com/webinars

*“One of the key features of organoids in general and those used in this study is the defined media that we use. This allows us to tailor the *in vitro* niche environment to the specific cells we are growing or to the cells that we are trying to generate in culture. This has been made possible, in part, through the use of a wide range of different growth factors and small molecules from R&D Systems and Tocris. We also used Cultrex RGF BME for all of the organoids in this study.”*

Dr. Talya Dayton, Hubrecht Institute, The Netherlands.

INTESTINAL ORGANOIDS

The small intestine, large intestine, and colon consist of a multicellular epithelium with distinct morphological structures, including villi and invaginated crypt structures. Intestinal crypts house Lgr5+ intestinal adult stem cells that are responsible for the continuous renewal of intestinal epithelium and were first utilized to create long-term 3D culture models of the intestine, termed intestinal organoids or epithelial organoids. These organoid cultures are employed to study normal and diseased physiology, including barrier functions, nutrient uptake, and tissue renewal. In addition, intestinal organoids can be generated from iPSCs. iPSC-derived organoids have been used as advanced models for gastrointestinal developmental biology, drug toxicity, and personalized medicine applications.

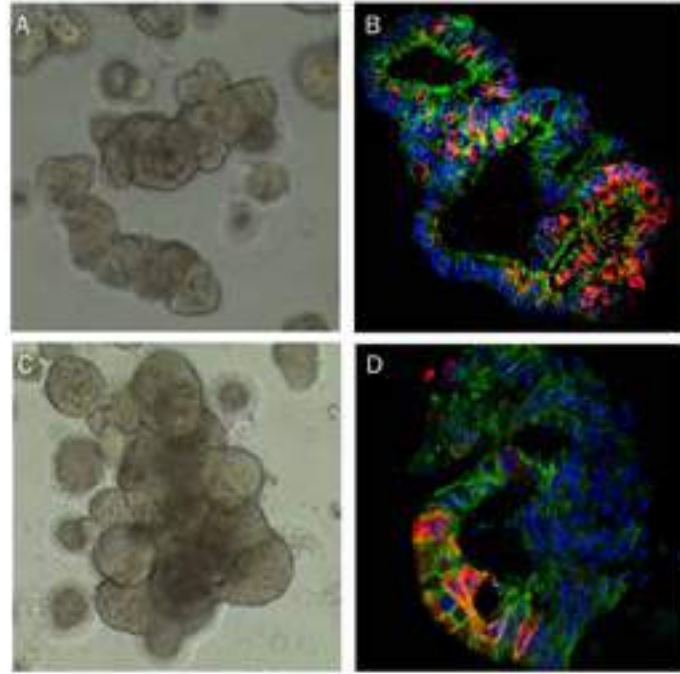

Human Intestinal Organoids Cultured using Cultrex RGF BME, Type 2. Human transverse colon organoids (A, B) and human ileum organoids (C, D) were grown using cells isolated from transverse colon and ileum biopsy tissue, respectively. Organoids were embedded in **Cultrex RGF BME, Type 2** (R&D Systems, Catalog # 3533-005-02) as a scaffold matrix. **A)** Brightfield image of human transverse organoids. **B)** Human transverse organoid stained using a **Goat Anti-Human/Mouse E-Cadherin Antigen Affinity-purified Polyclonal Antibody** (green; R&D Systems, Catalog # AF748), a **Mouse Anti-Human MUC2 Monoclonal Antibody** (red; Novus Biologicals, Catalog # NBP2-44431), and **DAPI** (blue; R&D Systems, Catalog # 5748). **C)** Brightfield image of human ileum organoids. **D)** Human ileum organoid stained using a **Rabbit Anti-Human Aldolase B Polyclonal Antibody** (red; Novus Biologicals, Catalog # NBP2-15345), a **Mouse Anti-Human Cadherin-17 Monoclonal Antibody** (green; R&D Systems, Catalog # MAB1032), and **DAPI** (blue; R&D Systems, Catalog # 5748).

NOTABLE PUBLICATIONS AND PROTOCOLS

PUBLICATION	DESCRIPTION	BIO-TECHNE REAGENTS USED	CATALOG #
Mahe M. <i>et al.</i> (2013) <i>Curr. Protoc. Mouse Biol.</i> 3 :217.	Establishment of gastrointestinal epithelial organoids.	Recombinant Mouse Noggin	6997-NG
		Recombinant Mouse Wnt-3a	1324-WN
		Recombinant Human EGF	236-EG
		Recombinant Human Jagged 1 Fc Chimera Protein	1277-JG
Co, J.Y. <i>et al.</i> (2019) <i>Cell Reports</i> 26 :2509	Human enteroid model for host-pathogen interactions.	Cultrex RGF BME, Type 2	3533-005-02
		A 83-01	2939
		CHIR 99021	4423
Sato T. and H. Clevers (2015) <i>Cell</i> . 161 :1700	Review of protocols for growing organoids from stem cells.		
Jung P. <i>et al.</i> (2011) <i>Nature Medicine</i> . 17 :1225	Isolation of human colonic stem cells.		

REAGENTS USED FOR INTESTINAL ORGANOID CULTURE

PRODUCT NAME	CATALOG #
Cultrex Reduced Growth Factor Basement Membrane Extract (RGF BME), Type 2	3533-005-02
GlutaminePlus	B90210
HEPES	R35150
N21-MAX Supplement	AR008
N-2 MAX Supplement	AR009
N-Acetylcysteine	5619
Gastrin I (Human)	3006
Y-27632 dihydrochloride (Rho Kinase Inhibitor)	1254
Recombinant Human EGF	236-EG
Recombinant Human R-Spondin 1	4645-RS
Recombinant Human Noggin	6057-NG
Recombinant Human FGF-2	3718-FB
A 83-01	2939
CHIR 99021	4423
Recombinant Human Wnt-3a	5036-WN

"Our research is greatly facilitated by Bio-Techne products. The various organoid systems and co-cultures are all performed in or on R&D Systems Cultrex Reduced Growth Factor BME with great results."

Jens Puschof, Hans Clevers lab, Hubrecht Institute, The Netherlands.

GASTRIC ORGANOIDS

Similar to the intestine, the stomach contains Lgr5+ adult stem cells that can be isolated, cultured, and differentiated *in vitro* into gastric organoids. Early organoid models elucidated molecular mechanism underlying gastric development, including signaling pathways that influence fundic or antral gastric epithelium formation. Gastric organoid cultures are powerful models to study normal and diseased gastric physiology as well as more complex models for drug discovery and disease modeling.

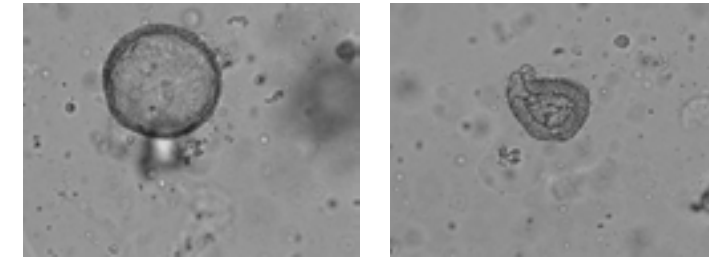


FIGURE 1. Undifferentiated Human Gastric Organoids. Representative brightfield images of human gastric organoids that were cultured using Cultrex RGF BME, Type 2 and the Bio-Techne reagents listed in this protocol.

REAGENTS USED FOR GASTRIC ORGANOID CULTURE

REAGENT NAME	SUPPLIER	CATALOG #
Cultrex Organoid Harvesting Solution	R&D Systems	3700-100-01
Cultrex Reduced Growth Factor Basement Membrane Extract (RGF BME) Type 2	R&D Systems	3533-005-02
Glutamine	Tocris Bioscience	5823
HEPES	Tocris Bioscience	3173
N21-MAX Supplement	R&D Systems	AR008
N-2 MAX Supplement	R&D Systems	AR009
N-Acetylcysteine	Tocris Bioscience	5619
Gastrin I (Human)	Tocris Bioscience	3006
SB 202190 (p38 MAPK Inhibitor)	Tocris Bioscience	1264
Nicotinamide	Tocris Bioscience	4106
Human Insulin, Solution	Sigma-Aldrich	I9278
Human Transferrin	Sigma-Aldrich	T8158
Y-27632 dihydrochloride (Rho Kinase Inhibitor)	Tocris Bioscience	1254
Recombinant Human EGF	R&D Systems	236-EG
Recombinant Human R-Spondin 1	R&D Systems	4645-RS
Recombinant Human Noggin	R&D Systems	6057-NG
Recombinant Human FGF-10	R&D Systems	345-FG
A 83-01 (ALK5 Inhibitor)	Tocris Bioscience	2939
CHIR 99021, (GSK-3 Inhibitor)	Tocris Bioscience	4423
Recombinant Human Wnt-3a	R&D Systems	5036-WN

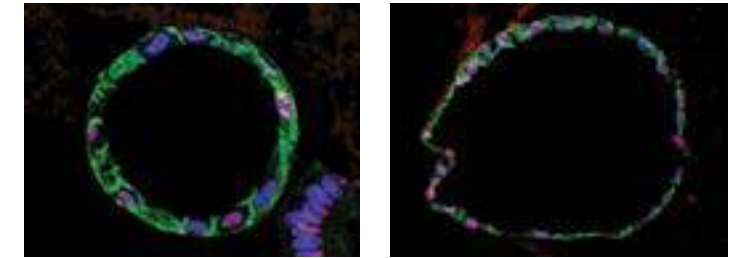



FIGURE 2. Immunohistochemistry of Undifferentiated Human Gastric Organoids. Human gastric organoids were cultured using Cultrex RGF BME, Type 2 and Bio-Techne reagents listed in this protocol. Undifferentiated colon organoids were stained using the Human/Mouse E-Cadherin Antibody (green; R&D Systems; Catalog # AF748), the Human HOXB7 Antibody (red; R&D Systems; Catalog # MAB8040), and counterstained with DAPI (blue; Tocris; Catalog # 5748).

NOTABLE PUBLICATIONS AND PROTOCOLS

PUBLICATION	DESCRIPTION	BIO-TECHNE REAGENTS USED	CATALOG #
McCracken, E.M. <i>et al.</i> (2014) <i>Nature</i> 516 :400. McCracken, E.M. <i>et al.</i> (2017) <i>Nature</i> 541 :182.	Modelling human development and disease in pluripotent stem-cell-derived gastric organoids.	Recombinant Human BMP4	314-BP
		Recombinant Human FGF4	235-F4
		Recombinant Human Noggin	6057-NG
		Recombinant Human EGF	236-EG
		Recombinant Human FGF10	345-FG
		Recombinant Human Wnt5a	645-WN
Munero J. <i>et al.</i> (2017) <i>Cell Stem Cell</i> 21 :51.	iPSC differentiation into colonic organoids.	Y-27632 dihydrochloride	1254
		Recombinant Human FGF4	235-F4
		Recombinant Human BMP2	355-BM
		Recombinant Human EGF	236-EG
		Recombinant Human Noggin	6057-NG
		CHIR 99021	4423
		SAG	4366
Li, X. <i>et al.</i> (2018) <i>Nature Comm.</i> 9 :2983.	Organoid cultures provide a model for recapitulate esophageal adenocarcinoma.	Cultrex BME RGF Type 2	3533-005-02
		Cultrex HA-R-Spondin1-Fc 293T Cells	3710-001-01
		A83-01	2939

LIVER ORGANOIDS

The liver is the primary organ system for drug metabolism and detoxification. In this role, it is also highly susceptible to damage from pharmaceuticals and other chemical toxicants. Animal models and traditional *in vitro* assays modeling liver metabolism often fail to recapitulate the *in vivo* toxicity of drugs in human patients. Liver organoids, derived from primary tissue or induced pluripotent stem cells, have emerged as more complex and predictive models for hepatotoxicity and drug screening.

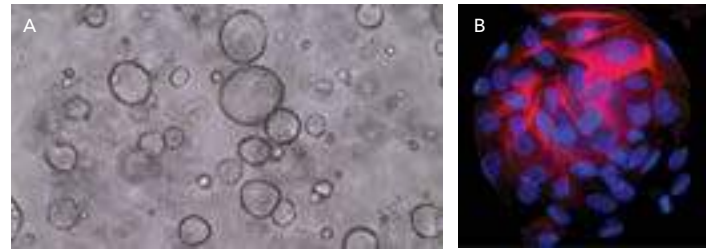

Human Liver Organoids. A) Brightfield Image of human undifferentiated liver organoids cultured using Cultrex RGF BME, Type 2 and in media featuring Bio-Techne reagents. B) Expression of Albumin (red; R&D Systems, Catalog # MAB1455) in differentiated human liver organoids. Image counterstained with DAPI (blue; Tocris, Catalog # 5748).

REAGENTS USED FOR LIVER ORGANOID CULTURE


REAGENT NAME	SUPPLIER	CATALOG #
Cultrex Organoid Harvesting Solution	R&D Systems	3700-100-01
Cultrex Reduced Growth Factor Basement Membrane Extract (RGF BME), Type 2	R&D Systems	3533-005-02
Glutamine	Tocris Bioscience	5823
HEPES	Tocris Bioscience	3173
N21-MAX Supplement	R&D Systems	AR008
N-2 MAX Supplement	R&D Systems	AR009
N-Acetylcysteine	Tocris Bioscience	5619
Gastrin I (Human)	Tocris Bioscience	3006
Nicotinamide	Tocris Bioscience	4106
Y-27632 dihydrochloride (Rho Kinase Inhibitor)	Tocris Bioscience	1254
Recombinant Human EGF	R&D Systems	236-EG
Recombinant Human R-Spondin 1	R&D Systems	4645-RS
Recombinant Human Noggin	R&D Systems	6057-NG
Recombinant Human FGF-10	R&D Systems	345-FG
Recombinant Human FGF-19	R&D Systems	969-FG
Recombinant Human BMP7	R&D Systems	354-BP
Recombinant Human HGF	R&D Systems	294-HG
Forskolin	Tocris Bioscience	1099
A 83-01 (ALK5 inhibitor)	Tocris Bioscience	2939
Recombinant Human Wnt-3a	R&D Systems	5036-WN
DAPT	Tocris Bioscience	2634
Dexamethasone	Tocris Bioscience	1126

NOTABLE PUBLICATIONS AND PROTOCOLS

PUBLICATION	DESCRIPTION	BIO-TECHNE REAGENTS USED	CATALOG #
Huch, M. <i>et al.</i> (2015) Cell 160:299.	Long-term culture of adult human liver stem cells.	Cultrex RGF BME, Type 2	3533-010-02
		A 83-01	2939
		Forskolin	1099
Ogawa, M. <i>et al.</i> (2015) Nat. Biotechnol. 33:853. cc	Human iPSC-derived cholangiocyte organoids.	L-685,458	2627
		Recombinant Human HGF	294-HG
		Recombinant Human EGF	236-EG
		Recombinant Human TGF-β1	240-B
Takebe, T. <i>et al.</i> (2013) Nature 499:481.	Vascularized human iPSC-derived hepatic organoids.	Recombinant Human Oncostatin M (OSM)	295-OM
Broutier, L. <i>et al.</i> (2016) Nature Protocols 11:1724.	Protocol for generating human and mouse adult liver organoids.	Cultrex RGF BME, Type 2	3533-005-02
Koike, H. <i>et al.</i> (2019) Nature 574:112.	Modelling hepato-biliary-pancreatic organogenesis.		

LUNG ORGANOIDS

3D cell culture models of the pulmonary system are increasingly utilized to study lung regeneration, model disease (*i.e.* cystic fibrosis), and investigate mechanisms of viral lung infection (*i.e.* SARS-CoV-2). While lung organoids were first generated using Lgr5⁺ stem cells isolated from primary tissue, protocols for culturing iPSC-derived lung organoids have increased the flexibility and accessibility of this model system for use in personalized medicine and drug discovery.

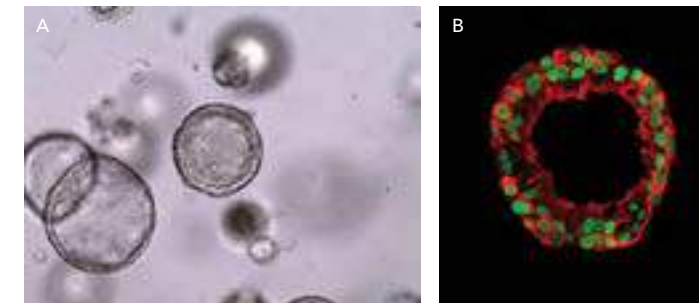

Human Lung Organoids. A) Brightfield Image of human lung organoids cultured using Cultrex RGF BME, Type 2 and in media featuring Bio-Techne reagents. B) Expression of Sox2 (green; R&D Systems, Catalog # AF2018) and Acetylated Tubulin (red; Novus Biologicals, Catalog # NB600-567) in human lung organoids.

REAGENTS USED FOR LUNG ORGANOID CULTURE


REAGENT NAME	SUPPLIER	CATALOG #
A 83-01	Tocris Bioscience	2939
Cultrex Organoid Harvesting Solution	R&D Systems	3700-100-01
Cultrex Reduced Growth Factor Basement Membrane Extract (RGF BME), Type 2	R&D Systems	3533-005-02
Advanced DMEM/F-12 Cell Culture Medium		
Glutamine	Tocris Bioscience	5823
HEPES	Tocris Bioscience	3173
N21-MAX Supplement	R&D Systems	AR008
N-Acetylcysteine	Tocris Bioscience	5619
Penicillin/Streptomycin	R&D Systems	B21210
SB 202190 (p38 MAPK Inhibitor)	Tocris Bioscience	1264
Nicotinamide	Tocris Bioscience	4106
Y-27632 dihydrochloride (Rho Kinase Inhibitor)	Tocris Bioscience	1254
Recombinant Human R-Spondin 1	R&D Systems	4645-RS
Recombinant Human Noggin	R&D Systems	6057-NG
Recombinant Human FGF-10	R&D Systems	345-FG
Recombinant Human FGF-7	R&D Systems	251-KG

NOTABLE PUBLICATIONS AND PROTOCOLS

PUBLICATION	DESCRIPTION	BIO-TECHNE REAGENTS USED	CATALOG #
Sachs, N. <i>et al.</i> (2019) EMBO J. 34:e100300.	Protocol for adult stem cell-derived lung organoids.	Cultrex RGF BME, Type 2	3533-005-02
Miller, A. J. <i>et al.</i> (2019) Nature Protocols 14:518.	Protocol for generating human iPSC-derived lung organoids.	Recombinant Human Activin A	338-AC
		Recombinant Human Noggin	6057-NG
		Recombinant Human FGF-10	345-FG
		Recombinant Human FGF-4	7460-F4
		Recombinant Human FGF-7	251-KG
Dye, B.R. <i>et al.</i> (2015) eLife 4:e05098.	Protocol for generating human iPSC-derived lung organoids.	Recombinant Human Activin A	338-AC
		Recombinant Human Noggin	6057-NG
		Recombinant Human FGF2	233-FB
		Recombinant Human FGF4	7460-F4
		Recombinant Human Sonic Hedgehog	8908-SH

BRAIN ORGANOIDS

Protocols to generate 3D brain organoids from ESCs and iPSCs were first published in 2009. These studies showed that pluripotent stem cells could differentiate into cerebral organoids containing specific cortical regions, neural progenitor populations, and cortical layer patterning. Cerebral organoids have since been employed to uncover evolutionary differences in brain development between species, mechanisms of brain region interconnectivity, and the developmental physiology of normal and diseased brain regions. iPSC-derived organoids show great potential for use in drug discovery as well as modeling neurodegenerative disease and viral brain infection.

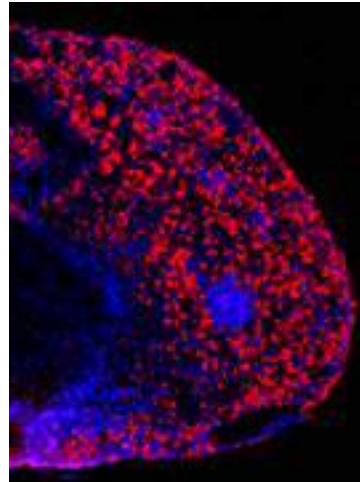

Neural Progenitor Cells in Cerebral Organoids. Human iPSCs were differentiated into cerebral organoids following the Lancaster protocol. Organoids were harvested and stained using a Human Pax6 Polyclonal Antibody (Catalog # AF8150) to identify neural progenitor cells in the developing cerebral tissue. Tissue was counterstained using DAPI (Catalog # 5748).

NOTABLE PUBLICATIONS AND PROTOCOLS


PUBLICATION	DESCRIPTION	BIO-TECHNE REAGENTS USED	CATALOG #
Pollen, A.A. <i>et al.</i> (2019) <i>Cell</i> 176 :743.	Establishing Cerebral Organoids as Models of Human-Specific Brain Evolution	Y-27632 dihydrochloride	1254
		SB 431542	1614
Bershteyn, M. <i>et al.</i> (2017) <i>Cell Stem Cell</i> 20 :435.	Human iPSC-Derived Cerebral Organoids Model Cellular Features of Lissencephaly and Reveal Prolonged Mitosis of Outer Radial Glia	Y-27632 dihydrochloride	1254
		SB 431542	1614
Bagley, J. A. <i>et al.</i> (2017) <i>Nature Methods</i> 13 :743.	Fused cerebral organoids model interactions between brain regions		
Lancaster M.A. and J. A. Knoblich (2014) <i>Nat. Protocols</i> 9 :2329.	Generation of cerebral organoids from human pluripotent stem cells		

REAGENTS USED FOR BRAIN ORGANOID CULTURE

PRODUCT NAME	CATALOG #
BASE MEDIA COMPONENTS	
N-2 MAX Supplement	AR009
N21-MAX Supplement	AR008
N21-MAX Vitamin A Free Supplement	AR012
Penicillin/Streptomycin	B21210
GlutaMAX	B90210
Insulin	
2-mercaptoethanol	
MATRIX AND MEDIA ADDITIVES	
Cultrex RGF Basement Membrane Extract, Type 2	3533-005-02
Recombinant Human FGF basic	3718-FB
Recombinant Human Noggin	6057-NG
Y-27632 dihydrochloride	1254

KIDNEY ORGANOIDS

Using pluripotent stem cells, kidney organoid culturing protocols have shown the ability to recapitulate the organ's complex tissue cytoarchitecture, including expression of cellular markers for podocytes, proximal tubules, and distal tubules. Success in cultivating kidney organoids has facilitated research interrogating kidney development, physiology, and mechanisms underlying kidney disease (*i.e.* chronic kidney disease). In addition, kidney organoid research has demonstrated its potential as a translational method for kidney tissue regeneration.

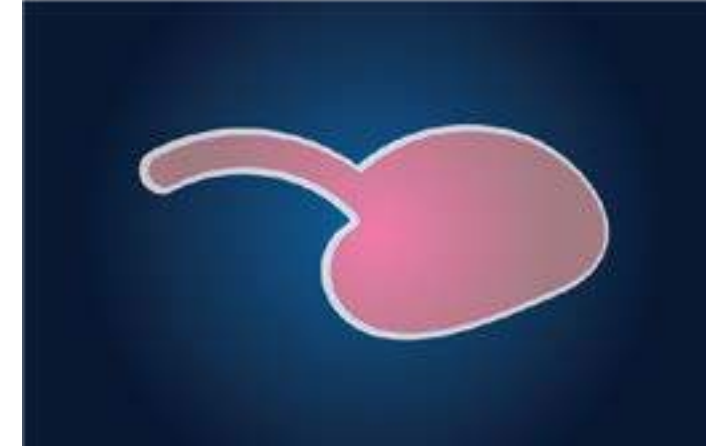

REAGENTS USED FOR KIDNEY ORGANOID CULTURE

PRODUCT NAME	CATALOG #
BASE MEDIA COMPONENTS	
N-2 MAX Supplement	AR009
N21-MAX Supplement	AR008
N-Acetylcysteine	5619
Penicillin/Streptomycin	B21210
GlutaMAX	B90010
Holo-Transferrin	
Advanced DMEM/F-12	
MATRIX AND MEDIA ADDITIVES	
Cultrex RGF Basement Membrane Extract, Type 2	3533-005-02
Recombinant Human Activin A	338-AC
Recombinant Human BMP-2	355-BM
Recombinant Human BMP-4	314-BP
Recombinant Human FGF basic	3718-FB
Recombinant Human FGF-9	273-F9
CHIR 99021	4423
Retinoic Acid	695
Y-27632 dihydrochloride	1254

NOTABLE PUBLICATIONS AND PROTOCOLS

PUBLICATION	DESCRIPTION	BIO-TECHNE REAGENTS USED	CATALOG #
Taguchi, A. <i>et al.</i> (2014) <i>Cell Stem Cell</i> 14 :53.	iPSC-nephron progenitor role in developing kidney organoids.	Recombinant Human Activin A	338-AC
		Recombinant Human BMP-4	314-BP
		Recombinant Human FGF-9	273-F9
		Recombinant Human FGF-2	233-FB
Morizane, R. <i>et al.</i> (2015) <i>Nat. Biotechnol.</i> 33 :1193.	Nephron organoids derived from human pluripotent stem cells.	Recombinant Human Activin A	338-AC
		Recombinant Human FGF-9	273-F9
		Y-27632	1254
		CHIR 99021	4423
Freedman, B.S. <i>et al.</i> (2015) <i>Nat. Comm.</i> 6 :8715.	Gene editing of kidney organoids to model disease.	IWP 2	3533
Takasato, M. <i>et al.</i> (2015) <i>Nature</i> 526 :564.	Human iPSC-derived kidney organoid generation.		

HEART ORGANOIDS

In vitro generation of cardiac tissue is enabling advancements in drug discovery and toxicity testing, as well as facilitating the engineering of cardiac tissue for regenerative therapies. Various methods have been employed to generate 3D cardiac tissue, including iPSC-derived cardiomyocyte spheroids and bioprinting of cardiac organoids with iPSCs that are subsequently differentiated into cardiomyocytes. However, protocol and reagent advancements are still needed to enhance the maturity and complexity of the cardiac tissue.



NOTABLE PUBLICATIONS AND PROTOCOLS

PUBLICATION	DESCRIPTION	BIO-TECHNE REAGENTS USED	CATALOG #
Andersen, P. (2018) Nat. Comm. 9:3140.	Precardiac spheroids generated from human pluripotent stem cells.	Recombinant Human Activin A	338-AC
		Recombinant Human BMP-4	314-BP
		Recombinant Human Wnt-3a	5036-WN
		Recombinant Human Wnt-5a	645-WN
		Recombinant Human Wnt-11	6179-WN
Kupfer, M.E. et al. (2020) Circulation Research 127:207.	Cardiac organoid formation using different extracellular matrix proteins.		
Mills, R.J. et al. (2017) PNAS 113:E8372.	Cardiac organoids from human iPSCs		

MAMMARY ORGANOIDS

Protocols to generate mammary organoids from primary epithelial tissues are helping elucidate the cell fate decisions and molecular mechanisms of mammary gland development, including ductal formation and transformation of milk-producing alveoli. Most importantly, these 3D culture techniques have enabled the cultivation of a breast cancer organoids, which are being employed for *in vitro* drug discovery and personalized drug screening for breast cancer.




NOTABLE PUBLICATIONS AND PROTOCOLS

PUBLICATION	DESCRIPTION	BIO-TECHNE REAGENTS USED	CATALOG #
Sachs, N. et al. (2018) Cell 172:373.	Robust protocol for long-term culturing of human mammary epithelial organoids.	Cultrex RGF BME, Type 2	3533-005-02
		Recombinant Human R-Spondin 3	3500-RS
Rosenbluth, J.M. (2020) Nat. Comm. 1711.	Human mammary organoids derived from breast tissue.	A83-01	2939
Jamieson, P. et al. (2017) Stem Cells and Regeneration 144:1065.	Mouse mammary organoids derived from epithelial cells.	Cultrex RGF BME, Type 2	3533-005-02
		Y-27632 dihydrochloride	1254

PANCREATIC ORGANOIDS

Pancreatic organoids have become an informative *in vitro* model to study pancreatic cancer, exocrine disease, and the basic development of pancreatic ductal epithelium for potential use as regenerative or therapeutic treatment of diabetes. While robust protocols for pancreatic organoid generation using mouse primary pancreatic ductal tissues exist, protocols that support the long-term cultivation of pancreatic organoids from human tissues are still emerging.

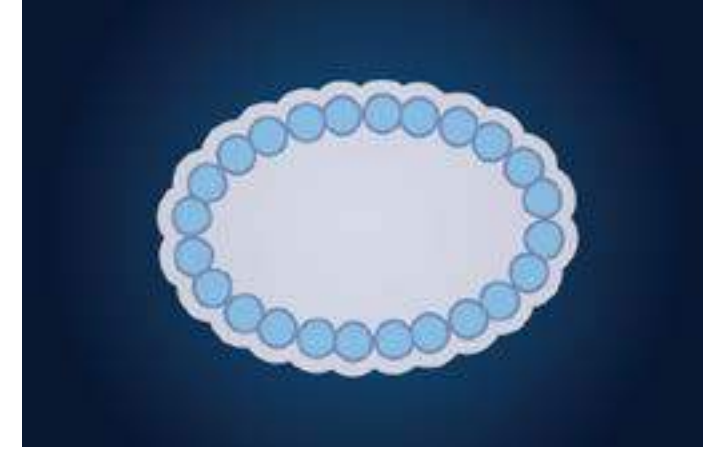

REAGENTS USED FOR PANCREATIC ORGANOID CULTURE


PRODUCT NAME	CATALOG #
BASE MEDIA COMPONENTS	
N-2 MAX Supplement	AR009
N21-MAX Supplement	AR008
N-Acetylcysteine	5619
Penicillin/Streptomycin	B21210
Glutamine	B90010
Advanced DMEM/F-12	
MATRIX AND MEDIA ADDITIVES	
Cultrex RGF Basement Membrane Extract, Type 2	3533-005-02
Recombinant Human EGF	236-EG
Recombinant Human FGF-10	345-FG
Recombinant Human Noggin	6057-NG
Recombinant Human R-Spondin 1	4645-RS
Recombinant Human Wnt-3a	5036-WN
A 83-01	2939
Nicotinamide	4106
Gastrin	3006

NOTABLE PUBLICATIONS AND PROTOCOLS

PUBLICATION	DESCRIPTION	BIO-TECHNE REAGENTS USED	CATALOG #
Georgakopoulos, N. et al. (2020) BMC Developmental Biology 20:4.	Long-term expansion of adult human pancreatic organoids.	Cultrex RGF BME, Type 2	3533-005-02
		PGE-2	2296
		Forskolin	1099
		A83-01	2939
Broutier, L. et al. (2016) Nat. Protocols 11:1724.	Protocol to culture self-renewing human pancreatic organoids.	Cultrex RGF BME, Type 2	3533-005-02
		Cultrex HA-R-Spondin1-Fc 293T Cells	3710-001-01
		Recombinant Human FGF-19	2939
		A83-01	969-FG
Greggio, C. et al. (2013) Development 140:4452.	Effects of reagents and matrices on pancreatic organoid culture.	R-spondin 1	4645-RS
		FGF-10	345-FG
		FGF-2	233-FB
Dossena, M. et al. (2020) Stem Cell Research & Therapy 11:94.	GMP-compliant culture of human pancreatic organoids.	R-spondin 1	4645-RS
		A83-01	969-FG
		PGE-2	2296

INNER EAR ORGANOIDS

Pluripotent stem cell-derived inner ear organoids are rapidly advancing our understanding of inner ear development and physiology. Inner ear organoids have been shown to develop sensory epithelium containing the necessary hair cells, supporting cells, and synaptic-like structures that support auditory or gravitational transduction. These models have great potential for translational research, uncovering molecular and cellular mechanisms that support the regeneration of cochlear and vestibular sensory tissue.

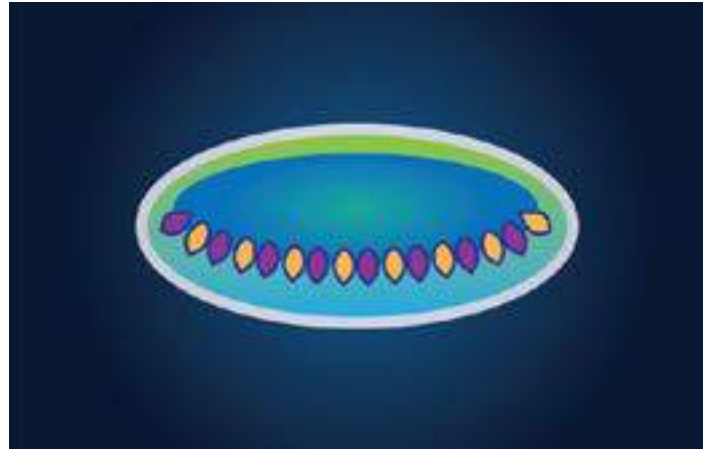

NOTABLE PUBLICATIONS AND PROTOCOLS

PUBLICATION	DESCRIPTION
Koehler K.R. and E. Hashino (2014) Nat. Protocols. 9 :1229.	Inner ear organoids from mouse embryonic stem cells.
Liu, X. <i>et al.</i> (2016) Nat. Comm. 7 :11508.	Vestibular organoids from mouse embryonic stem cells.
Koehler, K.R. <i>et al.</i> (2017) Nat. Biotech. 35 :583.	Inner ear organoids from human induced pluripotent stem cells.

REAGENTS USED FOR INNER EAR ORGANOID CULTURE

PRODUCT NAME	CATALOG #
LIF-2: BASE MEDIA COMPONENTS	
N-2 MAX Supplement	AR009
N21-MAX Supplement	AR008
Leukemia Inhibitor Factor (LIF)	7734-LF
CHIR 99021	4423
Penicillin/Streptomycin	B21210
Glutamine	B90010
Advanced DMEM/F-12	
MATRIX AND MEDIA ADDITIVES	
Cultrex RGF Basement Membrane Extract, Type 2	3533-005-02
Recombinant Human BMP-4	314-BP
Recombinant Human FGF basic	3718-FB
A 83-01	2939

ORGANOID HARVESTING

Organoids are often cultured in matrix hydrogels that promote the growth of 3D structures but also must be removed before passaging, cryopreservation, and analysis of the organoids. Proteases can be employed to degrade the extracellular proteins within the organoid matrix. However, non-enzymatic methods of matrix depolymerization, such as Cultrex Organoid Harvesting Solution, are preferred because they limit carryover of protease activity in subsequent cultures or product analysis.

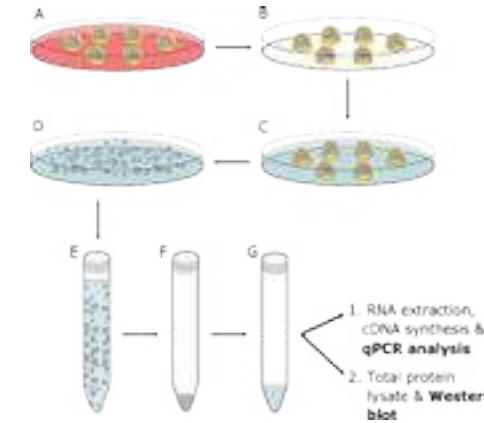


Figure 1. Summarized Protocol to Harvest Organoids for Biochemical Analysis. A) Treat organoids with differentiation medium. B) Discard medium. C) Add Cultrex Organoid Harvesting Solution. D) Incubate at 2-8 °C. E) Transfer organoids to a conical tube. F) Centrifuge organoids. G) Resuspend organoids in the appropriate lysis solution for either RNA extraction or protein analysis.

PRODUCT NAME	CATALOG #
Cultrex Organoid Harvesting Solution	3700-100-01

ORGANOID CRYOPRESERVATION

Cryopreservation of organoids is useful for cell line banking or when generating repositories of patient-derived organoids for drug discovery or toxicology testing. Similar techniques and reagents used to freeze down cell lines and primary cells can be employed for organoid cryopreservation, including base medium containing 20% FBS and 10% DMSO*. Due to their complex structural elements, troubleshooting cell viability during cryopreservation is a technical challenge. Freezing media and freeze-down strategy may need to be customized by tissue-type, organoid maturation, structure (freezing of intact structures, partially dissociated fragments, or as fully dissociated single cell suspensions), and density. R&D Systems offers basic reagents currently being used to prepare organoid cryopreservation media.

PRODUCT NAME	CATALOG #	BRAND	DESCRIPTION
Fetal Bovine Serum - Premium Select	S11150	R&D Systems	USDA-approved, Central American-origin FBS
Fetal Bovine Serum - Optima	S12450	R&D Systems	US-Origin, USDA-APHIS approved FBS
DMSO, sterile filtered	3176	Tocris	

*Reference organoid cryopreservation protocol based on Taniguchi Laboratory at MD Anderson Cancer Center.



WHERE SCIENCE
INTERSECTS INNOVATION™

bio·techne®

bio-techne.com

R&D SYSTEMS

**NOVUS
BIOLOGICALS**

TOCRIS

protein
simple

A&D

exosome_x

Global info@bio-techne.com bio-techne.com/find-us/distributors TEL +1 612 379 2956 North America TEL 800 343 7475
Europe | Middle East | Africa TEL +44 (0)1235 529449 China info.cn@bio-techne.com TEL +86 (21) 52380373

For research use or manufacturing purposes only. Trademarks and registered trademarks are the property of their respective owners.

BR_Organoid Handbook_STRY0049616